Tips: Before use, use a small amount of water to wipe the conductive sheet or wet towel to wipe the use site!



**USER MANUAL** 

### Instructions for use

#### Product model: JT-66

Product name: Cervical Massager

Rated power: 3W Charging time is about 2~3 hours Charging parameters: 5V === 1A

Rated voltage: 3.7V Rated time: 15min

Executive standard: GB4706.1-2005 GB4706.10-2008

#### JT-66 Charging operation instructions


- 1. 

  power button, long press power on/off to turn on/off the massage function.

  2. Power button, short press the adjustment mode.
- There are 6 modes, each time you switch, the corresponding gear will be emitted and the corresponding
- sound will be emitted.

  3. ⊕ pulse intensity upshift, a total of 15 levels.
- 3. ⊕ pulse intensity upshift, a total of 15 levels.
  4. There are 15 levels of ⊝ -pulse force reduction, and a beep will sound every time you switch.
  5. Charging function, use 5V = 1A mobile phone charging head to charge, the charging red light will be on, and it will go out when it is full.
  6. Heating control. The heating is turned on by default when the machine is turned on. Short press the heating button to turn off/on the heating. button to turn off/on the heating.

# Conductive Sheet/ heating Sheet

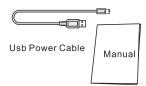

## People who should not use this product

(1) Take good care of children to ensure that they do not play with this

product.
(2) The product has a hot surface, so people who are not sensitive to heat

(2) The product has a hot surface, so people who are not sensitive to heat must pay attention to it when using it.

(3) Patients with allergic constitution heart disease and other people with electronic devices in their bodies should be used with caution.

(4) This product is not intended to be used by people (including children) who are physically impaired, sensory or mentally impaired, or lack relevant experience and knowledge, unless someone responsible for their safety supervises or directs their use of this product.

If you encounter problems when using this product, please refer to the following tips and possible solutions.

Status	Possible cause/treatment		
Not working	●Please confirm that the switch is turned on.  ●Please check whether the battery has power.		
Suddenly the massage stopped	Please check whether the product battery is exhausted.     Inadvertently press the button on the remote control.     The 10 minutes set by the automatic timing device has passed.     If the problem still exists, please do not continue to use and contact the company's maintenance service center.		
The machine shuts down automatically	When there is no contact with the skin, the machine will make a sound and it will automatically shut down after 1 minute. When the battery is low, the machine will automatically shut down after 10 beeps. Please charge		

# Safety need to know

When using this massage product, especially when children are present, the following safety precautions must be observed. To avoid the risk of electric shock from rabbits, please read all instructions carefully before use.

\*\*Be sure to pull out the power plug after use and before wiping.

\*\*Do not use when taking a shower or bathing. If this appliance falls into water, please pull out the plug immediately.

\*\*Do not place this product in a place where it can easily slip or fall into the bath and sink.

\*\*Do not use in a humid environment, do not put this product in water or other liquids.

WDo not use needles or other metal objects to fix this product, do not use unrecommended accessories and accessories, and do not use it in the case that the appliance shell or power cord is damaged.

WIf the product is accidentally wetted into the water, please do not use it by

#If the product is accidentally wetted into the water, please do not use it by yourself.

#Do not use this appliance when you are drunk or feel unwell. Do not use this product within one hour after a meal.

#If you have a neck injury or have recently undergone surgery, please use it under the guidance of a doctor.

#This product is used under battery power supply. If the battery is not charged or the battery is damaged, stop using it immediately, send it to the designated maintenance outlet or send it back to the manufacturer for battery replacement or maintenance. The user must not replace the battery without permission!

#Before discarding this product, take out the battery as follows: knock the product casing with a hammer along the edge, find the location where the battery is placed, and then use scissors to cut off the battery leads one by one and take out the battery (be careful when using scissors);

#Before discarding the product, the battery must be removed from the product;

product;

—When removing the battery, the product must be disconnected from the grid power supply. The battery should be safely disposed of.

If the product is still in a fault state or cannot be used, please turn off the power and unplug it and contact our service center.

Users should not disassemble or repair this product by themselves.

Touching the internal parts may cause malfunctions or electric shock accidents. Product damage caused by self-disassembly or repair will invalidate the product warranty, and the company is not responsible for any consequences caused by this situation.

#### Warranty card

	Customer Name: Detailed Address:					
Customer Information	Correspondence address:		Postal code:			
Intormation	Contact phone: Mobile phone:					
Shopping Mall Name				Purchase Price		

#### Certificate of conformity

Product number: Inspector code: Check 03

Production Date: